

SenSura[®] Mio Baby

Long wear time - less disturbance

SenSura[®] Mio Baby

SenSura Mio Baby
delivers long wear time
without compromising
skin-friendliness.



There's a Mio for every body



The first SenSura® Mio pouching system designed specifically for premature babies and infants

SenSura® Mio Baby is a two-piece pouching system featuring skin friendly adhesive, which fits individual body profiles and covers multiple stomas. The specialized adhesive delivers long wear time and a secure fit to the body, while maintaining gentle application and removal from the skin.

SenSura Mio Baby is designed to ease your daily practice and:


- Reduce the time spent handling the infant
- Minimize disturbance for the infant
- Make it easier to meet individual needs (no starter hole)
- Provide access to the stoma

Additional infant-friendly features




Skin friendly adhesive

The specialized adhesive sticks well and delivers long wear time, protects skin against aggressive output that causes skin irritation, and removes gently using the Brava® Baby Adhesive Remover.

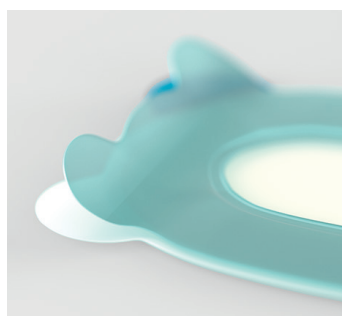


True 2-piece system

The adhesive coupling's soft flange is designed to be comfortable for babies, and easy to handle for caregivers - it's easy to apply and remove the pouch from the flange, and the oval shape enables the coupling to accommodate multiple stomas if needed.

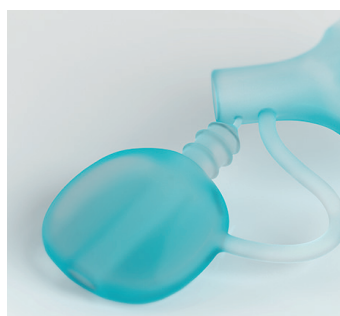

Soft, easy-to-use outlet

The outlet is compatible with enteral syringes for easy emptying while the soft material is gentle on baby's skin.

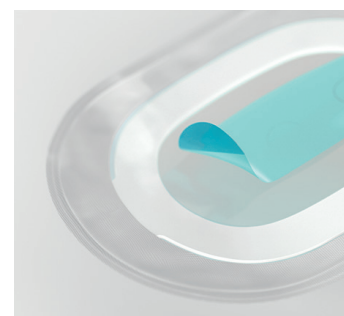

Access window

A dedicated access window facilitates easy access to the stoma and to refeed, for less disturbance of the baby.

Product assortment list

Item code	Product name	Coupling size	Units / box	HCPCS*
18700	2-piece Flex drainable pouch, with soft spout outlet and access window, neutral gray		10	A5063
18705	2-piece Flex barrier, oval cut-to-fit, 0 - 1 9/16"		5	A4414
12014	Brava® Baby Adhesive Remover (3 mL liquid sachet)	–	30	A4455

Prior to use, refer to the product 'Instructions for Use' for intended use and relevant safety information.

1. Compared to standard of care. Verified by laboratory testing.

*Reimbursement Disclaimer: Coloplast Corp. provides this information for your general reference and related to the reimbursement of Coloplast products only. Reimbursement, coverage and payment policies can vary from one insurer and region to another, and may change over time. Coloplast does not guarantee coverage or payment of products.

Ostomy Care / Continence Care / Wound & Skin Care / Interventional Urology

Coloplast Corp. Minneapolis, MN 55411 / 1-800-533-0464

www.coloplast.us The Coloplast logo is a registered trademark of Coloplast A/S. © 2022 Coloplast Corp. All rights reserved.



PM-05897 07.22